

Report

3355-28-0.qgd
D:\ding_ke\pql\3355-28-0.qgd
Admin 2015-9-21 10:54:42

Method

[Comment]

===== Analytical Line 1 =====

[GC-2010]

Column Oven Temp. :60.0 °C
Injection Temp. :220.00 °C
Injection Mode :Split
Flow Control Mode :Pressure
Pressure :93.1 kPa
Total Flow :79.5 mL/min
Column Flow :1.50 mL/min
Linear Velocity :44.7 cm/sec
Purge Flow :3.0 mL/min

Oven Temp. Program


Rate	Temperature(°C)	Hold Time(min)
-	60.0	0.00
10.00	100.0	4.00
10.00	200.0	20.00
10.00	250.0	10.00

\$EndIf\$

[GCMS-QP2010 Plus]

IonSourceTemp :250.00 °C
Interface Temp. :250.00 °C
Solvent Cut Time :1.00 min

Peak#	R.Time	I.Time	Peak Report TIC		Area	Area%	Height	
			F.Time				Height	Height%
1	2.252	2.217	2.483		11048074	100.00	3917567	100.00
					11048074	100.00	3917567	100.00

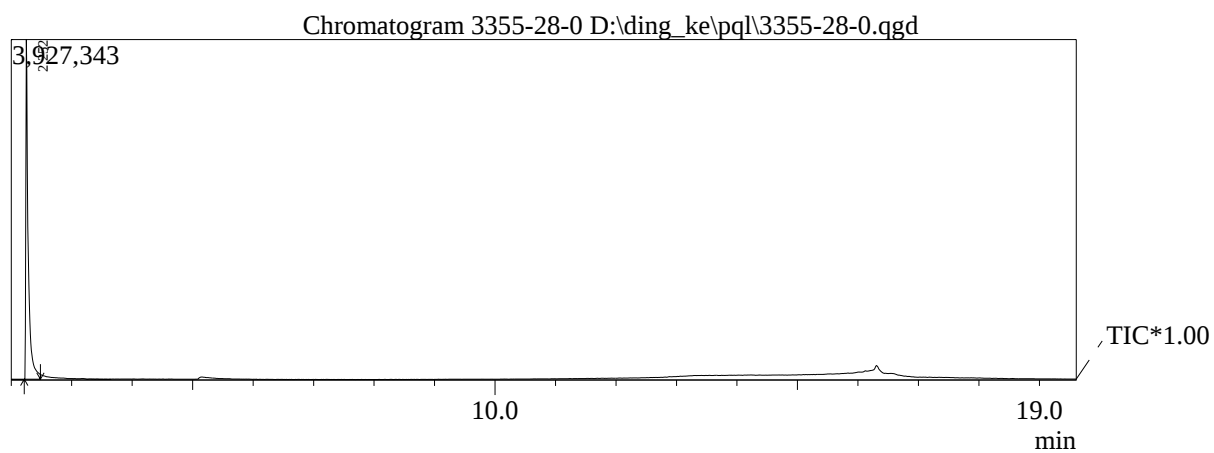

Spectrum

Line#:1 R.Time:2.3(Scan#:31)

MassPeaks:30

RawMode:Averaged 2.2-2.3(30-32) BasePeak:53(1725198)

BG Mode:Calc. from Peak Group 1 - Event 1

