

Report

90-02-8.qgd
D:\ding_ke\pql\90-02-8.qgd
Admin 2015-9-18 16:08:36

Method

[Comment]

===== Analytical Line 1 =====

[GC-2010]

Column Oven Temp. :60.0 °C
Injection Temp. :220.00 °C
Injection Mode :Split
Flow Control Mode :Pressure
Pressure :93.1 kPa
Total Flow :79.5 mL/min
Column Flow :1.50 mL/min
Linear Velocity :44.7 cm/sec
Purge Flow :3.0 mL/min

Oven Temp. Program


Rate	Temperature(°C)	Hold Time(min)
-	60.0	0.00
10.00	100.0	4.00
10.00	200.0	20.00
10.00	250.0	10.00

\$EndIf\$

[GCMS-QP2010 Plus]

IonSourceTemp :250.00 °C
Interface Temp. :250.00 °C
Solvent Cut Time :1.00 min

Peak#	R.Time	I.Time	Peak Report TIC F.Time	Area	Area%	Height	Height%
1	5.029	4.883	5.217	82354527	100.00	14297183	100.00
				82354527	100.00	14297183	100.00

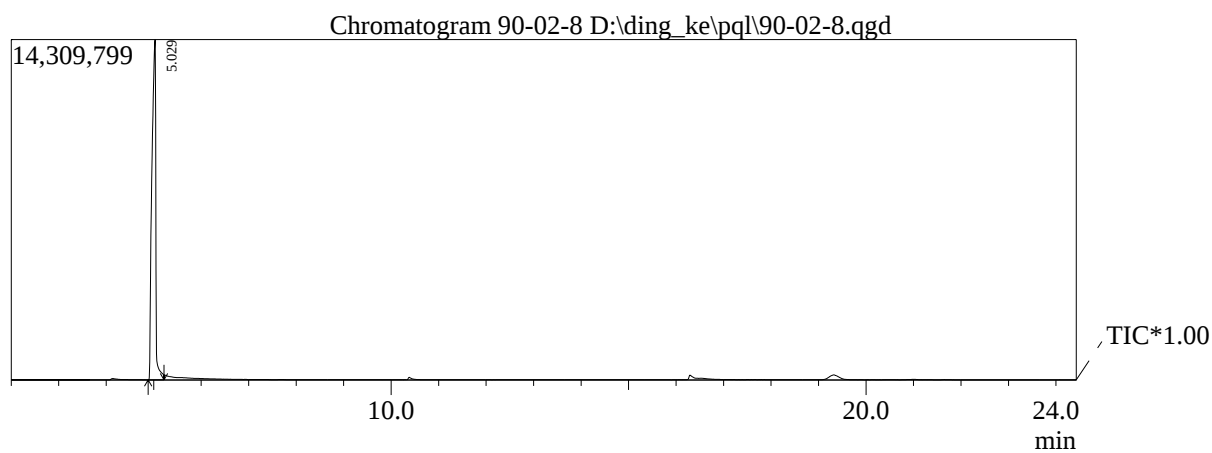

Spectrum

Line#:1 R.Time:5.0(Scan#:364)

MassPeaks:62

RawMode:Averaged 5.0-5.0(363-365) BasePeak:122(3678851)

BG Mode:Calc. from Peak Group 1 - Event 1

