

Report

5308-25-8.qgd
D:\ding_ke\pql\5308-25-8.qgd
Admin 2015-9-21 14:32:03

Method

[Comment]

===== Analytical Line 1 =====

[GC-2010]

Column Oven Temp. :60.0 °C
Injection Temp. :220.00 °C
Injection Mode :Split
Flow Control Mode :Pressure
Pressure :93.1 kPa
Total Flow :79.5 mL/min
Column Flow :1.50 mL/min
Linear Velocity :44.7 cm/sec
Purge Flow :3.0 mL/min

Oven Temp. Program


Rate	Temperature(°C)	Hold Time(min)
-	60.0	0.00
10.00	100.0	4.00
10.00	200.0	20.00
10.00	250.0	10.00

\$EndIf\$

[GCMS-QP2010 Plus]

IonSourceTemp :250.00 °C
Interface Temp. :250.00 °C
Solvent Cut Time :1.00 min

Peak#	R.Time	I.Time	Peak Report TIC F.Time	Area	Area%	Height	Height%
1	3.580	3.492	3.617	99770055	92.74	25187698	93.80
2	11.378	11.333	11.533	7815922	7.26	1664571	6.20
				107585977	100.00	26852269	100.00

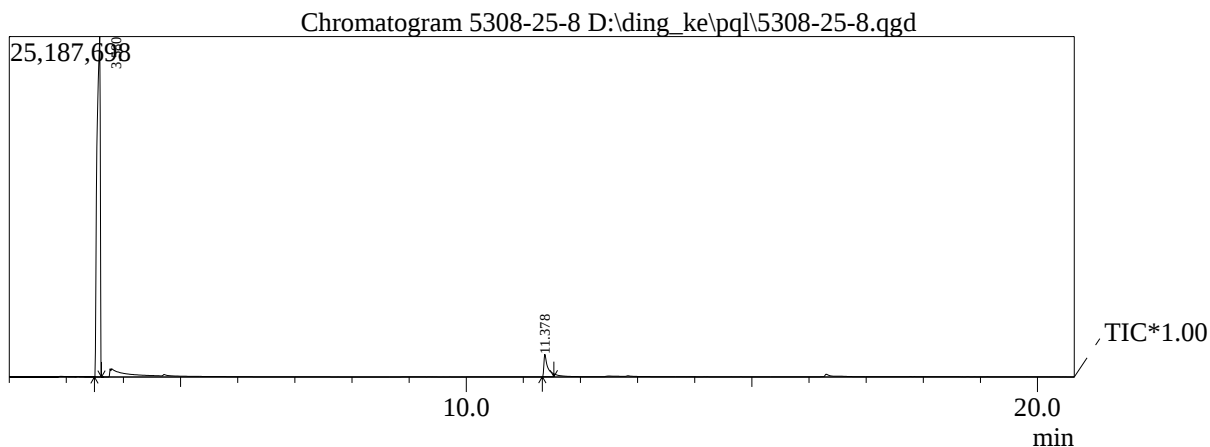


Spectrum

Line#:1 R.Time:3.6(Scan#:191)

MassPeaks:56

RawMode:Averaged 3.6-3.6(190-192) BasePeak:72(6323433)

BG Mode:Calc. from Peak Group 1 - Event 1

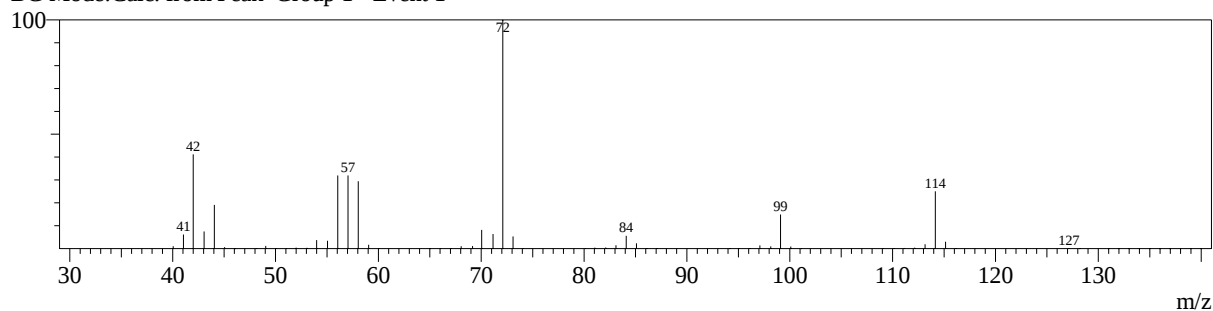


Line#:2 R.Time:3.6(Scan#:191)

MassPeaks:56

RawMode:Averaged 3.6-3.6(190-192) BasePeak:72(6323433)

BG Mode:Calc. from Peak Group 1 - Event 1

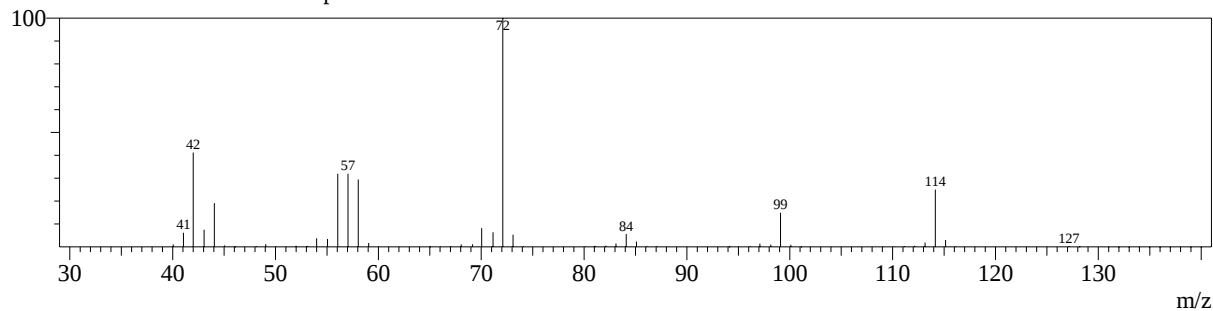

Line#:3 R.Time:11.4(Scan#:1126)

MassPeaks:49

RawMode:Averaged 11.4-11.4(1125-1127) BasePeak:127(172787)

BG Mode:Calc. from Peak Group 1 - Event 1

